CARRERA

Elektrische Zahnbürste

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Deutsch

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English

Manual
from page 12

Français

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Nederlands

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Elektrische Zahnbürste

Electric toothbrush

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Elektrische tandenborstel

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Dear customer,

We are pleased that you have chosen this product. If you have any questions, please contact customer service via our email address:

kundenservice@hsphanseshopping.de

Please read the instructions carefully before using the product for the first time and keep them in a safe place. If the product is passed on, these instructions must also be handed over. The manufacturer and importer accept no liability if the information in these instructions is not observed!

In the context of continuous further development, we reserve the right to change the product, packaging or instruction leaflets at any time.

SYMBOLES

Read carefully and follow the instructions to avoid personal injury and damage to property.

Read operating instructions before use!



Use indoors only!

The CE marking confirms that the appliance complies with the applicable European requirements that apply to the respective group of appliances and that the specific declarations of conformity are available.

SCOPE OF DELIVERY

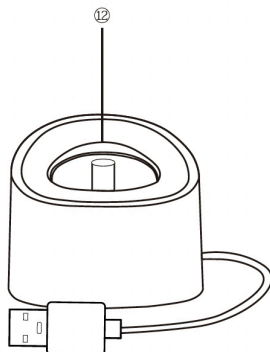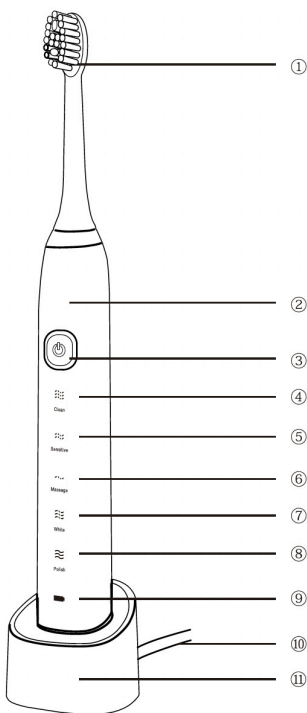1. Attachable brush
2. Handpiece
3. On / Off button
4. Cleaning mode „Clean“
5. Cleaning mode „Sensitive“
6. Cleaning mode „Massage“
7. Cleaning mode „White“
8. Cleaning mode „Polish“
9. Charging indicator
10. Power cable
11. Charging station
12. Handpiece holder

GENERAL NOTES

- Before using the unit for the first time: Carefully remove the unit from its packaging and check it for completeness and damage.
- The unit is for personal use only and is not intended for commercial applications.
- Only use the appliance as described in the instructions. Any other use is considered contrary to the intended use.
- The unit has an IPX7 protection class and is therefore protected against water jets from any angle and against the ingress of water for a limited time, but not for permanent submersion. The unit can be used in the bath without concern.

To ensure the maximum life of the product, clean it regularly as described in the instructions.

IMPORTANT SAFETY INSTRUCTIONS

- This appliance may also be used by children aged from 8 years and above and persons with reduced physical or sensory capabilities provided they have been given supervision and instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Ensure that children do not play with the appliance.
- Cleaning and maintenance by the user must not be carried out by children.
- CAUTION Danger of suffocation! Keep the packaging material away from children and animals.
- Check occasionally whether the unit/mains cable has any damaged spots. A damaged or defective unit/power supply unit must no longer be used. In this case, contact the customer service. Do not make any modifications to the unit. Otherwise there may be a risk of fire, electric shock or injury.
- Never pull the plug out of the socket by the cable. Do not touch the plug with wet hands. Otherwise there is a risk of electric shock.
- If you are undergoing dental treatment, you should contact your dentist before using this toothbrush.

- **WARNING:** There is a risk of explosion if the battery is not replaced correctly. Replace the battery only with the same or an equivalent type of battery.
- Always change the entire battery set. Do not mix old and new batteries or batteries with different charge levels.
- Never expose the battery to excessive heat (e.g. blazing sun) and never throw the battery into a fire.
- The connections must not be short-circuited.
- Keep the battery out of the reach of small children. If the battery is swallowed, seek medical attention immediately.
- Non-rechargeable batteries must not be recharged. Rechargeable batteries must be removed from the product before charging.
- If the battery has leaked, remove it from the battery compartment with the help of a cloth and dispose of it as intended. Avoid skin and eye contact with the battery acid. If you do have eye contact with the battery acid, rinse your eyes with plenty of water and contact a doctor immediately. In case of skin contact with the battery acid, wash the area with plenty of water and soap.
- Have damaged mains plugs and mains cables repaired or replaced immediately and only by authorised specialists!
- Never immerse the unit permanently in water or other liquids.
- Only use original accessories for this appliance.

OPERATION

Note: Replace the brushes every 3 months or sooner if they are worn out.

Commissioning

- Fully charge the unit before using it for the first time.
- To do this, connect the charging station to a designated socket using a commercially available USB adapter. Make sure that the charging station stands securely and cannot fall into the water.
- Place the handpiece on the handpiece holder.
- During the charging process, the charging indicator lights up red. When the charging indicator is green, the toothbrush is fully charged. If the charging indicator flashes red, the battery power is low and the appliance needs to be recharged.

Cleaning instructions

Moisten the brush head and apply any toothpaste of your choice. To avoid splashing the toothpaste, guide the brush to the teeth before switching on the appliance. Start with the outer sides first, then the inner sides and finally brush the chewing surfaces. Brush all four quadrants of the jaw evenly. To learn more about the correct brushing technique for you, ask your dentist or dental hygienist.

Note: The appliance has an automatic shut-off device which activates after two minutes and switches off the appliance.

Cleaning modes

Your toothbrush offers different cleaning modes that are displayed on the handpiece:

1. **CLEAN:** Standard cleaning mode for daily cleaning
2. **SENSITIVE:** Gentle and thorough cleaning for sensitive areas
3. **MASSAGE:** Gentle massage of the gums
4. **WHITE:** Brightening cleansing for occasional or daily use.
5. **POLISH:** Cleaning for a particularly clean mouthfeel (higher speed)

Press the corresponding button on the handpiece to switch between the cleaning modes.

CLEANING AND STORAGE

- Rinse the brush head thoroughly under running water for a few seconds after each use with the handpiece switched on.
- Always disconnect the mains plug before cleaning! Otherwise there is a risk of electric shock!
- Do not use abrasive, chemical or corrosive cleaning agents to clean the appliance.

TECHNICAL DATA

Device

Power supply: DC5V \equiv 1A

Protection class: IPX7/III

Battery

Typ: Li-Ion Akku 14500 3,7V 800mAh 2,96Wh

Charging time: approx. 10h

Operating time: up to 120 minutes

DISPOSAL

In connection with the batteries and rechargeable batteries placed on the market by us or with the delivery of devices containing batteries or rechargeable batteries, we are obliged to inform you of the following in accordance with the Battery Act (BattG): The symbol of the crossed-out dustbin means that batteries and accumulators must not be disposed of with household waste. If the battery or accumulator contains more than 0.0005% mercury by mass, more than 0.002% cadmium by mass or more than 0.004% lead by mass, the symbol outside the crossed-out dustbin is extended by the chemical symbols of the metals (Hg, Cd, Pb). As the end user, you are legally obliged not to dispose of all used batteries and accumulators in household waste, but to return them to retailers or municipal collection points. Used batteries and accumulators that are not enclosed in the old appliance, as well as lamps that can be removed from the old appliance, must normally be separated from the old appliance in a non-destructive manner before being handed in at a collection point. Old batteries and accumulators may contain harmful substances that can damage the environment or your health if not stored or disposed of properly. Returned old batteries and accumulators will only be recycled. This means that valuable raw materials such as iron, zinc, manganese or nickel are recovered. You can either return old batteries or accumulators free of charge in your immediate vicinity (e.g. to a retailer or recycling centre) or send them back to us free of charge. The environment and HSP Hanse Shopping GmbH thank you.

The symbol of the crossed-out dustbin means that this appliance is subject to the EU Directive on Waste Electrical and Electronic Equipment (WEEE Directive). It indicates that the respective appliance must be disposed of separately from unsorted municipal waste at the end of its life. A black bar under the crossed-out dustbin indicates that the appliance was placed on the market after 2005. This appliance may not be disposed of as normal household waste, but must be collected separately and disposed of in an environmentally sound manner at a municipal collection point. Owners of old appliances from private households can hand them in free of charge at the collection points of the public waste management authorities or at the take-back points set up by manufacturers or distributors within the meaning of the Electrical and Electronic Equipment Act. Shops with a sales area of at least 400 m² for electrical and electronic equipment are obliged to take back the equipment. From 01.07.2022, also those grocery shops with a total sales area of at least 800 m² which offer and make available on the market

electrical and electronic equipment several times a year or on a permanent basis are subject to this obligation. This also applies to sales using means of distance communication if the storage and dispatch areas for electrical and electronic equipment are at least 400 m² or the total storage and dispatch areas are at least 800 m².

The recycling symbol consists of three arrows that represent the reuse cycle (recycling). The packaging and filling materials can be recycled.

Do not dispose of the packaging and filling materials in household waste, but take them to the recycling collection in an environmentally friendly way.

WARRANTY

We assure that all products are made of high quality materials and have been checked with the utmost care. However, if damage or a fault should occur during normal use, please contact the seller of the product. A prerequisite for the warranty is operation and proper assembly in accordance with the operating instructions.

Improper use and / or improper transport may invalidate the warranty.

The warranty period is 2 years and begins on the date of purchase. If the product you have purchased is not free from defects, please contact the seller of the product within the period from the date of purchase.

The following are excluded from the warranty:

- Damage due to external force
- Non-compliance with the operating instructions
- Repairs and modifications by persons not authorised by us
- Damage due to incorrect handling or maintenance
- Damage due to neglect or accidents
- Wear and consumables

IMPORTANT NOTICE

If you have any questions about the product or warranty claims, please contact our customer service. Do not send your product to our address without request. The cost and risk of loss for unsolicited mailings will be borne by the sender. We reserve the right to refuse acceptance of unsolicited deliveries or to return corresponding goods to the sender at the sender's expense.

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