Sanelan[®]

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Vibrations-Zahnbürstenset 3tlg.

Artikel Nr. 85 33 72


Instruction Manual

3 Pcs Vibration Toothbrush Set

Article No. 85 33 72

CE

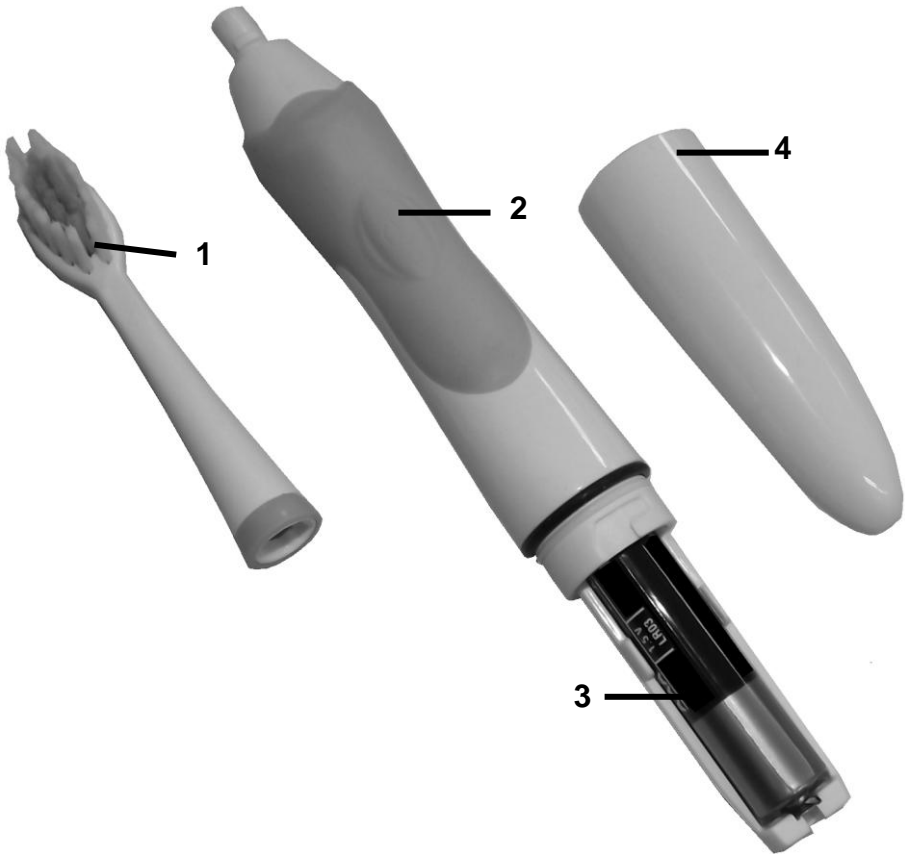1	Bürstenkopf / Brush Head
2	EIN/AUS-Taster / ON/OFF Button
3	Batteriefach /Battery Compartment
4	Batteriefachdeckel /Battery Compartment Lid

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
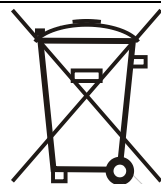
Safety Notes


Please note the following safety notes to avoid malfunctions, damage or physical injury:

- So not expose the toothbrush to high temperatures. The toothbrush is not a toy.
- Do not immerse the motor bloc in water.
- Do not disassemble the toothbrush and do not try to repair it yourself. Contact our customer support department.
- Do not disassemble, burn or short-circuit batteries.

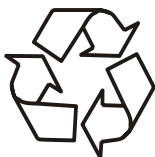Batteries do not belong in household garbage.

For battery disposal, please check with your local council.

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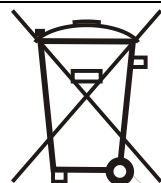Operation

Before first Use

Remove all parts from the packing. Dispose of packaging materials or store them out of reach of children.

Inserting the Battery

1. The battery compartment lid is locked with bayonet lock. Turn the battery compartment lid counter clockwise and remove it.
2. Insert one battery of the type LR03, 1.5 V (AAA) into the compartment. The battery symbol is marked inside the battery compartment.
3. Lock the battery compartment with the lid by turning clockwise.

Using

1. Attach the brush head to the motor unit. The brush head must face in the same direction as the ON/OFF button.
2. Before brushing your teeth, rinse your mouth thoroughly and wet the brush head. Use regular toothpaste.
3. To avoid splashing, switch on the toothbrush after guiding the brush head to your teeth. Place the brush in a slanted position.
4. Brush the rims of the gums, move the brush along the outside and the inside surfaces of all teeth, brush the front teeth and all molars. For further information

Operation

concerning the right cleaning technique, ask your doctor.

5. Brush your teeth with a rotary movement. You do not need to press hard to thoroughly remove plaque, and brush for at least 2 minutes. A change in the shape of the brush is normal after use. This is influenced by your individual cleaning technique, but does not affect the efficiency of the brush.
6. After you have finished brushing your teeth, switch off the unit. After cleaning your teeth, rinse your mouth thoroughly.

Cleaning

- Remove the brush head from the motor unit by turning counter clockwise.
- Rinse the brush head for several seconds under running water after use. Clean the motor unit with a towel. Never use aggressive chemicals or abrasives to clean the surfaces.
- Additionally, we recommend soaking the brush head in denture cleaner every 2 to 4 weeks. Replace the brush head after three months. Spare brush heads can be purchased from Westfalia under the article number 85 33 73.

Technical Data

Power Supply	1 Battery, Type LR03 (AAA)
Total Weight	40 g
Oscillation per Minute	22000
Degree of Protection	IPX4

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declare by our own responsibility that the product

3 Pcs Vibration Toothbrush Set

Article No. 85 33 72

is according to the basic requirements, which are defined in the European Directives

2004/108/EC

Electromagnetic Compatibility (EMC)

2011/65/EU

**Restriction of the Use of certain hazardous
Substances in electrical and electronic
Equipment (RoHS)**

and their amendments.

Hagen, March 10th, 2014

(Thomas Klingbeil, QA-Representative)

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Dear Customer,

Please help avoid waste materials.

If you at some point intend to dispose of this article, then please keep in mind that many of its components consist of valuable materials, which can be recycled.

Please do not discharge it in the rubbish bin, but check with your local council for recycling facilities in your area.